

ATI No. 164172  
ASTIA FILE COPY


STI

UNITED STATES ATOMIC ENERGY COMMISSION

NYO-2011


THE DETERMINATION OF NICKEL IN ZIRCONIUM  
WITH DIMETHYLGLYOXIME

By  
R. H. Beaumont, Jr.  
G. J. Petretic

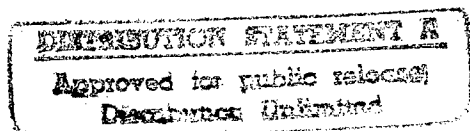

August 15, 1950

New Brunswick Laboratory

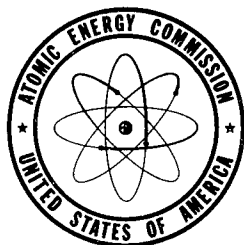

STANDARD INFORMATION

Technical Information Service, Oak Ridge, Tennessee

19961101 077

PRINTED IN USA  
PRICE 5 CENTS

Available from the  
Office of Technical Services  
Department of Commerce  
Washington 25, D.C.

## THE DETERMINATION OF NICKEL IN ZIRCONIUM WITH DIMETHYLGLYOXIME

By R. H. Beaumont, Jr. and G. J. Petretic

### ABSTRACT

A colorimetric method is presented for the determination of microgram quantities of nickel in zirconium metal with dimethylglyoxime. The zirconium metal is dissolved in dilute HF. Citric acid is added to complex the zirconium. Nickel is oxidized with bromine water. A large excess of ammonium hydroxide is added. The solution is cooled below room temperature and dimethylglyoxime is added. After approximately fifteen minutes, the optical density of the solution is measured at 540 m $\mu$ . Interferences are negligible. Results are given and are shown to be comparable to those obtained polarographically. The color is stable for at least three hours but results are better if the optical density is measured after a definite length of time.

\* \* \* \* \*

### I. Introduction.

A direct colorimetric method for the determination of nickel in zirconium metal has been requested. The method adopted is based on the reddish-brown color produced when tetravalent nickel reacts with dimethylglyoxime. The conditions are essentially those established by Pieters, Hanssen and Geurts (1) when determining nickel in iron. Slight modifications to further reduce errors due to interfering elements by using a narrow band as obtained by the Beckman spectrophotometer in place of a filter photometer have been made.

### II. Discussion.

Pieters et al (1) have succeeded in stabilizing the color of nickelic dimethylglyoxime, produced by bromine oxidation, by using a large excess of ammonium hydroxide. Consequently, the principal disadvantage of the method was overcome.

Pieters et al report that the principal interfering elements are Co, Fe, and Cr. In their work, however, they used a filter photometer with the center of the band width at approximately 520 m $\mu$ . When using a Beckman spectrophotometer with a 1 m $\mu$  band width at 538 m $\mu$ , the interference due to Fe and Cr is completely eliminated. They also report that Co interferes to

the extent of about 20% when present in a 10 to 1 ratio of Co to Ni. Since, however, only negligible quantities of cobalt have been found in the present zirconium production, the possibility of appreciable interference is neglected.

The sample is dissolved in dilute HF. Citric acid is added to complex the Zr, Al, Fe, etc., Bromine is used as the oxidizing agent. The solution is neutralized and an excess of ammonium hydroxide is added. The color is developed with dimethylglyoxime and its optical density is measured at 538 m $\mu$ , using a Beckman Spectrophotometer, 1 to 2 m $\mu$  band width and 1 cm cells.

### III. Experimental.

Pieters et al (1) report the use of concentrated NH<sub>4</sub>OH while in a separate part of the paper, NH<sub>4</sub>OH of density 0.95 is specified. The latter concentration corresponds to a more dilute solution and when used in the quantity suggested, a stable color is not established when determining nickel in zirconium. If, however, concentrated NH<sub>4</sub>OH is used in the same volume, a color which is nearly stable for at least three hours is produced. This is shown in Table 1.

Table 1.

#### STABILITY OF THE NICKELIC DIMETHYLGLYOXIME COMPLEX

0.5 g Zr/100 ml containing 5 ml excess conc. NH<sub>4</sub>OH. Measurements made at 540 m $\mu$  using 1 cm cells.

| $\gamma$ Added Ni | Optical Density |         |         |          |
|-------------------|-----------------|---------|---------|----------|
|                   | 15 min.         | 25 min. | 60 min. | 180 min. |
| 50                | .059            | .061    | .062    | .062     |
| 100               | .117            | .118    | .123    | .128     |
| 150               | .174            | .178    | .186    | .188     |
| 200               | .227            | .235    | .246    | .258     |
| 250               | -               | -       | .293    | .303     |

The method does not work well when more than 200  $\gamma$  ( $\approx$  400 ppm) of Ni are present. A precipitate forms and the nickelic dimethylglyoxime color complex is produced but slowly.

The color due to the ferric citrate complex is known not to interfere at wave lengths above about 500 m $\mu$  (2). The interference due to chromium as a chromate is also known to be negligible at 538 m $\mu$  (3). The possibility of chromium forming a color with dimethylglyoxime that might interfere at this wave length was investigated. No interference was detected due to chromium in varying amounts up to 500 ppm.

Inasmuch as the available standard zirconium samples (NBL 284 or NBL 285) were known to contain but little nickel, resort was made to other less positive means for determining the accuracy of the method. Samples were analyzed

for nickel which had previously been analyzed polarographically (4). A standard zirconium sample (NBL-285) was spiked with known quantities of nickel. Analysis of the spiked sample and three unknowns containing appreciable iron, chromium and nickel are shown in Table II.

Table II

COMPARISON OF COLORIMETRIC AND POLARO-  
GRAPHIC METHODS

0.500 gm. samples of Zirconium Metal

| <u>Sample No.</u> | <u>Ni Added-ppm</u> | <u>Ni Found - ppm</u> |                     |
|-------------------|---------------------|-----------------------|---------------------|
|                   |                     | <u>Polarograph</u>    | <u>Colorimetric</u> |
| 285               | 300                 | 290                   | 300                 |
| 6121              | 0                   | 65                    | 68                  |
| 6125              | 0                   | 140                   | 153                 |
| 6127              | 0                   | 310                   | 350                 |

#### IV. Procedure

1. Weigh a 0.500 gm sample of Zr metal into a 100 ml platinum or bakelite beaker.
2. Completely dissolve the sample in 10 ml of H<sub>2</sub>O plus 1 ml of 48% HF.
3. Add 10 ml of 10% citric acid, 5 ml of saturated Br<sub>2</sub> water, and about 5 ml of saturated boric acid (0.85M).
4. With stirring, add concentrated NH<sub>4</sub>OH until the color just disappears. Then add 5 ml in excess.
5. Transfer the solution to a 100 ml volumetric flask and cool to slightly below room temperature.
6. Add 3 ml of 0.1% dimethylglyoxime in alcohol.
7. Approximately 15 minutes after adding the dimethylglyoxime, dilute the solution to volume. Read the optical density of the solution in a 1.00 cm cell at 540 mμ.
8. Obtain the concentration in ppm by reference to a standard calibration curve.

#### SUMMARY

The dimethylglyoxime method for the determination of nickel has been applied to the analysis of zirconium samples. A stable color has been established after bromine oxidation. Results are accurate when determining amounts of nickel to about 200 γ.

## Bibliography

- (1) Pieters, Hanssen, and Geurts, Anal. Chim. Acta 2.377 (1948)
- (2) Sandell, "Colorimetric Determination of Traces of Metals", New York, Interscience Publishers, 1944.
- (3) Lingane and Collat, Anal. Chem., 22, 166 (1950)
- (4) Rodden, "Minutes of Zirconium Conference", NYO-2005, page 40, (1950)

END OF DOCUMENT